


TANT
Medical Engineering Ltd.

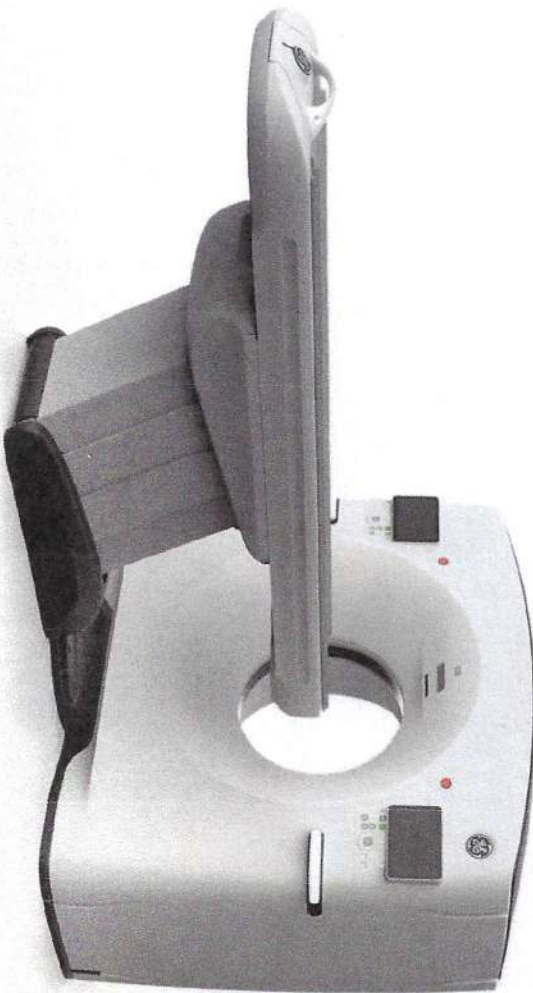

REV	A	24/FEB/2023	First issue drawing (DC-376449)
MODIFICATIONS			

- 01 - Cover Sheet
- 02 - Equipment Layout
- 03 - Structural - Electrical Layout
- 04 - Floor Structural Details
- 05 - Radiation Protection Layout
- 06 - Radiation Protection Details
- 07 - Power Requirements - Power Distribution
- 08 - Environment - Interconnections
- 09 - Equipment Dimensions (1)

- 10 - Equipment Dimensions (2)
- 11 - HVAC - Delivery
- 12 - Disclaimer - Site Readiness

SUDABELT MEDICAL CO LTD
LAGOS
NIGERIA



GE Healthcare
Efosa Amayo
Efosa.amayo@ge.com

REVOLUTION MAXIMA
FINAL STUDY

Drawn by	Verified by	Concession	S.O. (GON)	PIM Manual	Rev
N. Farfas	M. Czibók	-	-	5809942-1EN	6
Format	Scale	File Name		Date	Sheet
A3	1:50	CT-B329269-FIN-00-A.DWG		24/FEB/2023	01/12

A mandatory component of this drawing set is the GE Healthcare Pre Installation manual. Failure to reference the Pre Installation manual will result in incomplete documentation required for site design and preparation. Pre installation documents for GE Healthcare products can be accessed on the web at: www.gehealthcare.com/siteplanning

GE does not take responsibility for any damages resulting from changes on drawings made by others. Errors may occur by not referring to the complete set of final issue drawing. GE cannot accept responsibility for any damage due to the partial use of GE final issue drawings, however caused. All dimensions are in millimeters unless otherwise specified. Do not scale from printed pdf files. GE accepts no responsibility or liability for defective work due to scaling from these drawings.

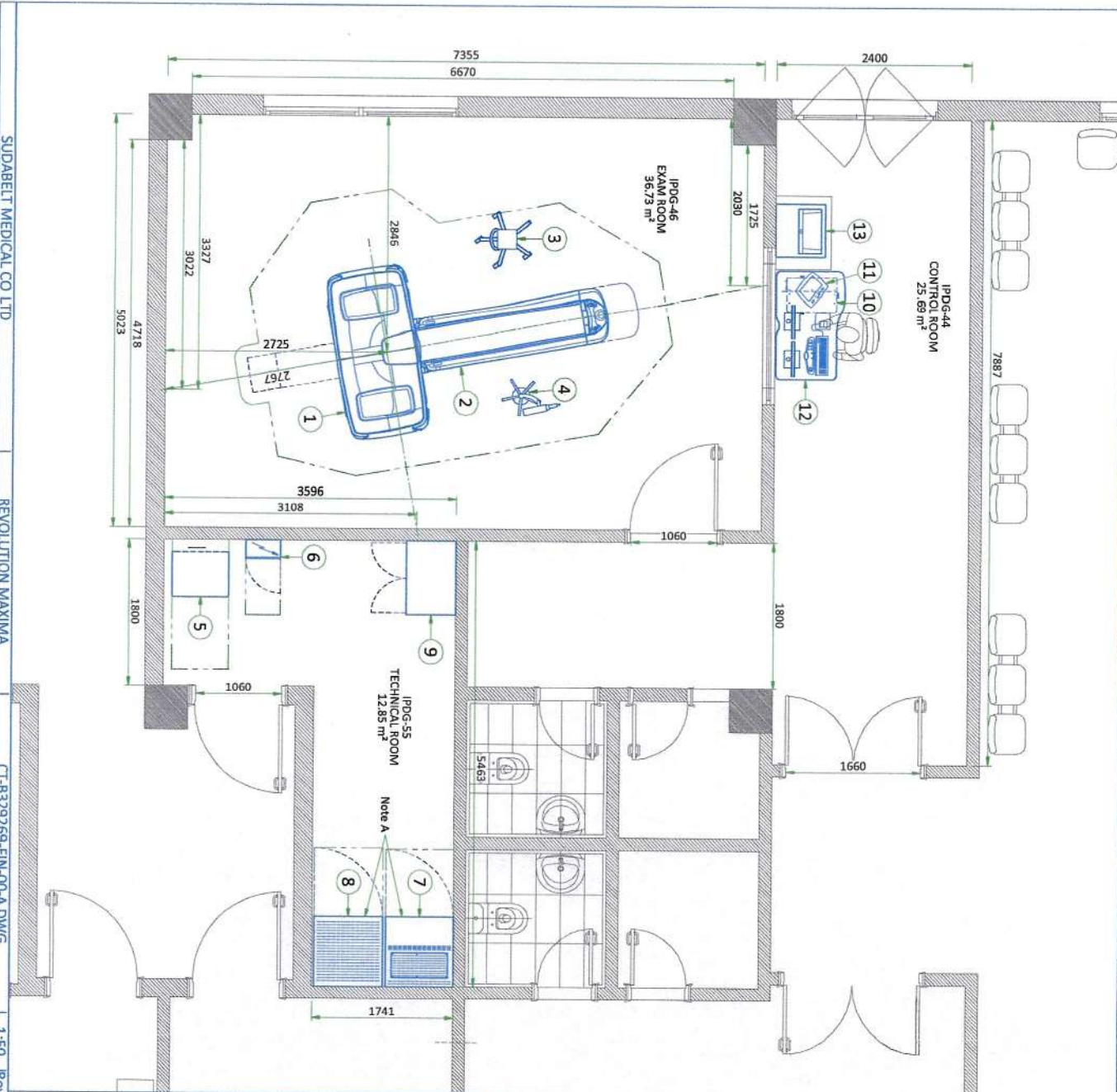

EQUIPMENT LAYOUT

ITEM	DESCRIPTION	DIMENSIONS LxWxH (mm)	WEIGHT (kg)
1	GANTRY	2041x1014x1928	1558
2	GT1700V TABLE	650x2370x1046	445
3	ECG MONITOR	-	3
4	INJECTOR ON PEDESTAL	-	-
5	POWER DISTRIBUTION UNIT (PDU)	700x550x1062	370
6	POWER DISTRIBUTION BOX (PDB)	225x424x929	33
7	FULL UPS 120 KVA	-	-
8	FULL UPS BATTERY CABINET	-	-
9	GE STORAGE CABINET	610x914x1067	68
10	OPEN CONSOLE	672x400x576	64.2
11	INJECTOR CONTROL	-	-
12	AURORA SWS TABLE	1300x850x850	40
13	DVS700 LASER IMAGER	660x610x470	54

WALL - ACCORDING TO RECEIVED DRAWING

STRUCTURE - ACCORDING TO RECEIVED DRAWING

EXAM ROOM HEIGHT


FINISHED FLOOR TO SLAB HEIGHT

FALSE CEILING HEIGHT

3.00 m
2.70 m

Notes:
A) Additional heat dissipation

Handwritten signature and date: 24/2/23

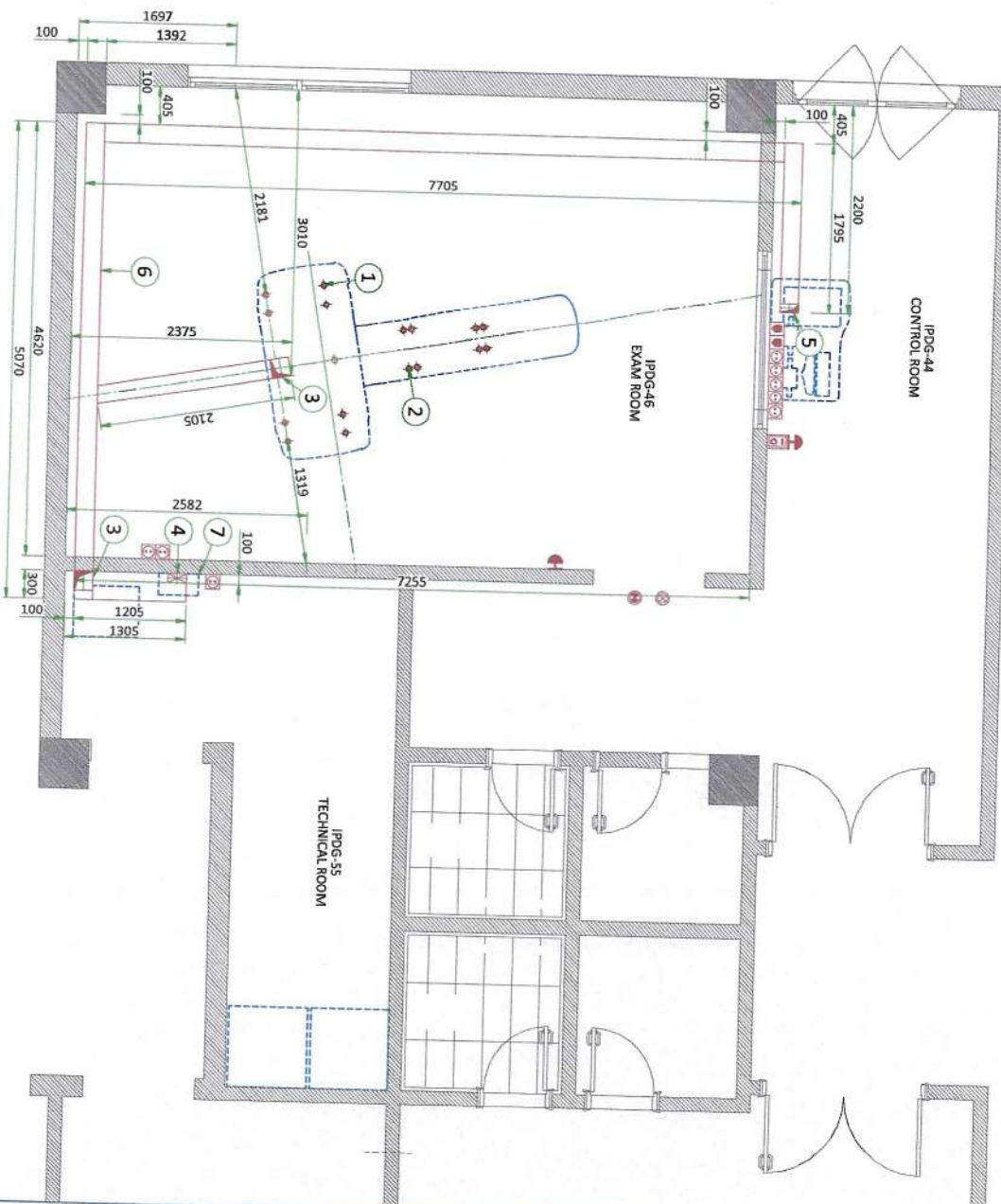


STRUCTURAL-ELECTRICAL LAYOUT

ITEM	QTY	DESCRIPTION
1		Gantry anchoring (see Structural Details)
2		Table anchoring (see Structural Details)
3		200x200 cable inlet on the floor
4		200x100 cable inlet on the floor and 200x100 vertical duct for MDP cabling (h=1.1m)
5		200x100 cable inlet on the floor
6		200x70 flush floor duct
7		Main Disconnect Panel (MDP)
		Electrical outlet 10/16A 230V + G
		FI 45 network socket
		System emergency off (SEO), (recommended height 1.50m-1.85m above floor)
		X-Ray ON lamp (L1) - 24V
		System ON light (L) - 24V
		System remote control (Y), locked when power OFF "ON" and "OFF" impulse buttons with indicator lamps red=ON / green=OFF located at 1.50m above floor
		Flush floor duct


Medical Engineering Ltd

ANCHORING/LOADING DISTRIBUTION TO THE FLOOR

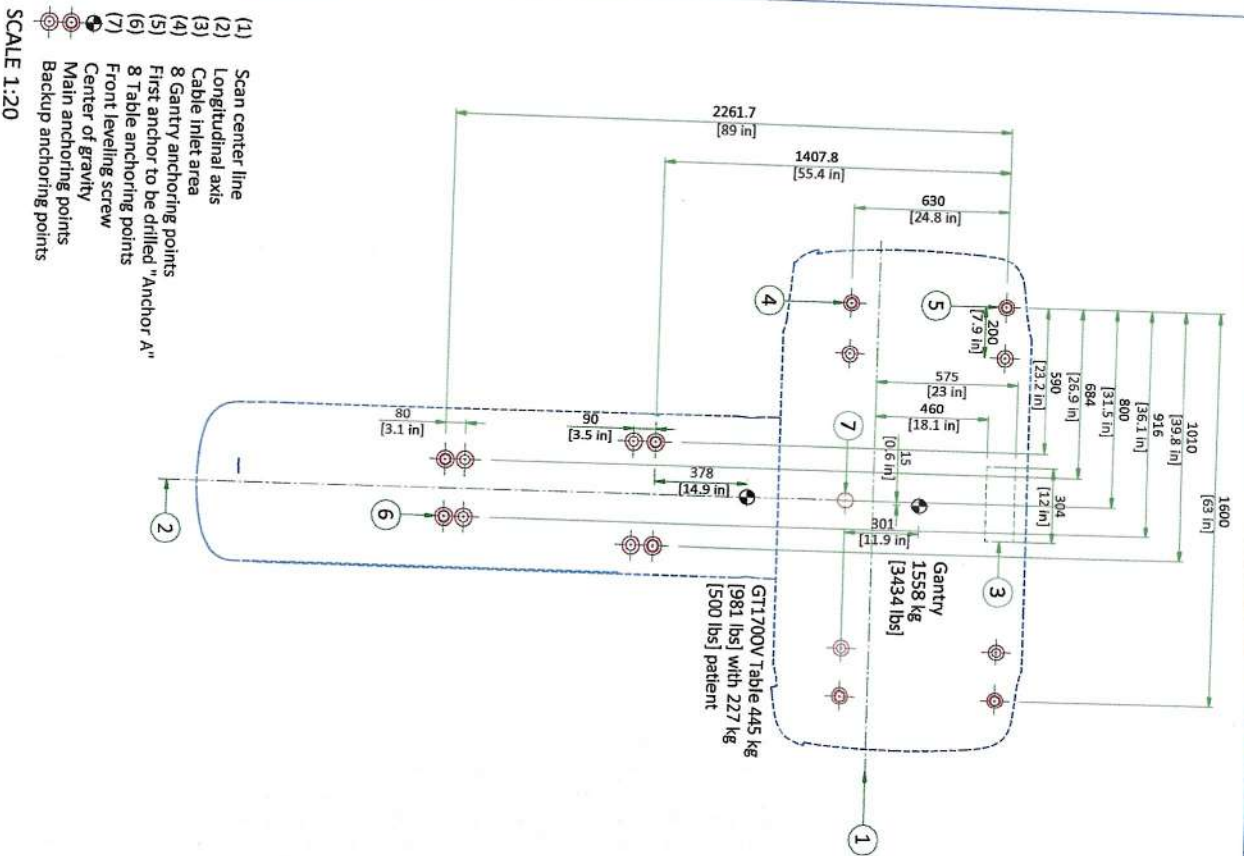

FLOOR REQUIREMENTS

GE SUPPLIED TABLE/GANTRY ANCHORS

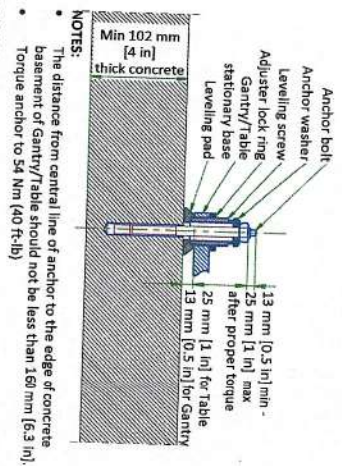


FINISHED FLOOR REQUIREMENTS

- Installation requires a finish floor in the scan and control rooms.
- The floor surface in the scan room directly under the gantry and table must be level.
- The floor levelness tolerance of the floor surface that the gantry and table will rest on is 6 mm [1/4 in] over a 3048 mm [10 ft] distance.
- Shims should not be used to compensate for a floor that does not meet this requirement.
- Eight or more floor covering openings that are 102 mm [4 in] in diameter are made to ensure the table and gantry rest on a solid surface.
- These requirements apply to all installation types.


NOT TO SCALE

TYPICAL CABLE MANAGEMENT

FLUSH FLOOR DUCT

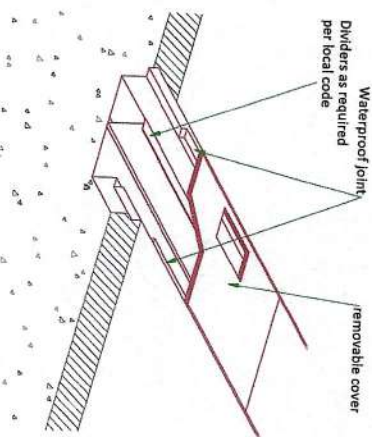

VERTICAL DUCT ON WALL

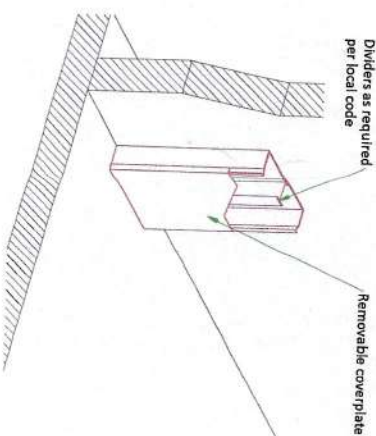

NOT TO SCALE

RADIATION PROTECTION LAYOUT

SHIELDING REQUIREMENTS SCALING

CHANGED PARAMETER (mas)	MULTIPLICATION FACTOR (new mas/100)
80 kV	0.24
100 kV	0.45
120 kV	0.71
140 kV	1.00
1 mm aperture	0.20
3 mm aperture	0.22
5 mm aperture	0.27
10 mm aperture	0.38
15 mm aperture	0.48
20 mm aperture	0.59
30 mm aperture	0.79
40 mm aperture	1.00

SHIELDING REQUIREMENTS:

Engage a qualified radiological health physicist to review your scan room shielding requirements, taking into consideration:


- Scatter radiation levels within the scanning room
- Equipment placement
- Weekly projected work-loads (number of patients/day technique (kVp*ma))
- Materials used for construction of walls, floor, ceiling, doors, and windows.
- Activities in surrounding scan room areas.
- Equipment in surrounding scan room areas (e.g., film developer, film storage)
- Room size and equipment placement within the room relative to room size.

The illustrations on this page depict measured radiation levels within the scanning room, while scanning a 32 cm CTDI phantom with the technique shown:

- 140 kV
- 100 mA
- 1 sec
- 40 mm

Use the mAs, kV and aperture scaling factors in the table shown here to adjust exposure levels to the scan technique used at the site.


NOTE: Actual measurements can vary. Expected deviations equals $\pm 15\%$, expect for the 5 mA and 1 mm techniques, where variations may be greater (up to a factor of 2), due to the inherent deviation in small values. The maximum deviation anticipated for tube output equals $\pm 40\%$.

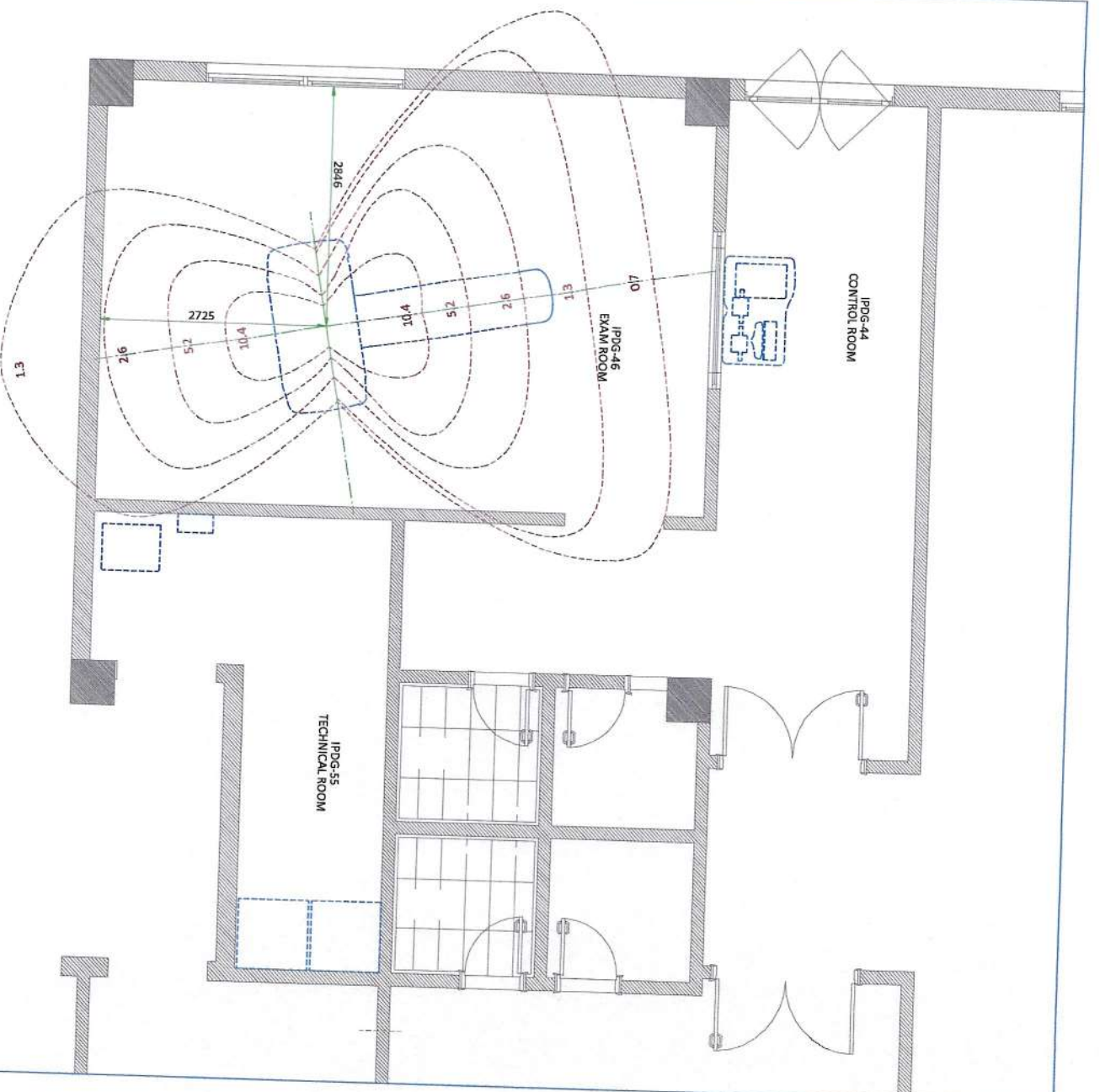

Michael J. McKeown, Ph.D.

RADIATION SCATTER - HEAD PHANTOM

NOTE: 140 kV
100 mAs/scan
1 sec
40mm aperture

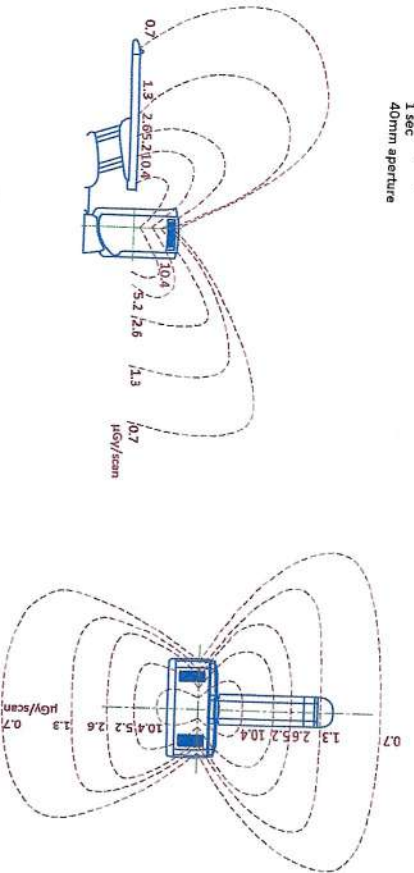


Elevation

Plan View

RADIATION SCATTER - BODY PHANTOM

NOTE: 140 kV
100 mAs/scan
1 sec
40mm aperture

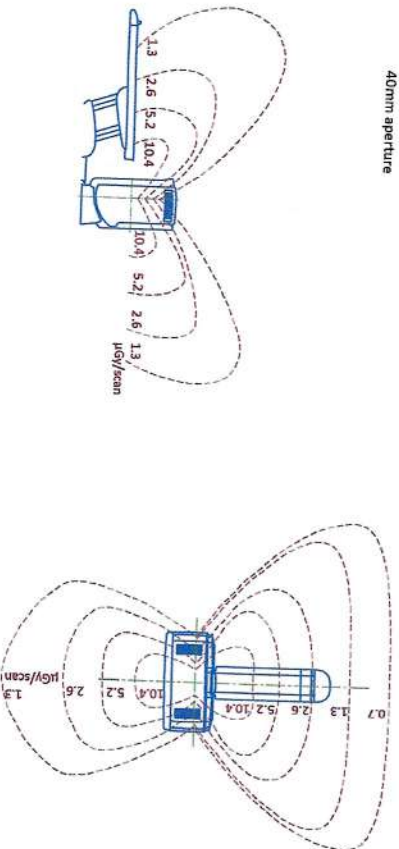

Elevation

Plan View

Handwritten signature and date: 24/2/2023

POWER REQUIREMENTS

POWER SUPPLY	3 PHASES+N+G
FREQUENCIES	200 /220/ 240/380/400/420/440/460/480 V ± 10%
MAXIMUM POWER DEMAND	50/60 Hz ± 3 Hz
AVERAGE (CONTINUOUS) POWER DEMAND	100 kVA
POWER FACTOR	20 kVA
	0.85

- Power supply should come into a main disconnect panel (MDP) containing the protective units and controls.
- The section of the supply cable should be calculated in accordance with its length and the maximum permissible voltage drops.
- There must be discrimination between supply cable protective device at the beginning of the installation (main low-voltage transformer side) and the protective devices in the MDP.

SUPPLY CHARACTERISTICS

- Power input must be separate from any others which may generate transients (elevators, air conditioning, radiology rooms equipped with high speed film changers...).
- All equipment (lighting, power outlets, etc...) installed with GE system components must be powered separately.
- Phase imbalance 2% maximum.
- Transients must be less than 1500V peak. (on a 400V line)

GROUND SYSTEM

- System of equipotential grounding.
- Equipotential: The equipotential link will be by means of an equipotential bar. This equipotential bar should be connected to the protective earth conductors in the ducts of the non GE cableways and to additional equipotential connections linking up all the conducting units in the rooms where GE system units are located.

CABLES

- Power and cable installation must comply with the distribution diagram.
- All cables must be isolated and flexible, cable color codes must comply with standards for electrical installation.
- The cables from signaling and remote control (V, SEO, L...) will go to MDP with a pigtail length of 1.5m, and will be connected during installation. Each conductor will be identified and isolated (screw connector).

CABLEWAYS

- The general rules for laying cableways should meet the conditions laid down in current standards and regulations, with regard to:
- Protecting cables against water (cableways should be waterproof).
- Protecting cables against abnormal temperatures (proximity to heating pipes or ducts).
- Replacing cables (cableways should be large enough for cables to be replaced).
- Metal cableways should be grounded.

Handwritten signature/initials

POWER DISTRIBUTION

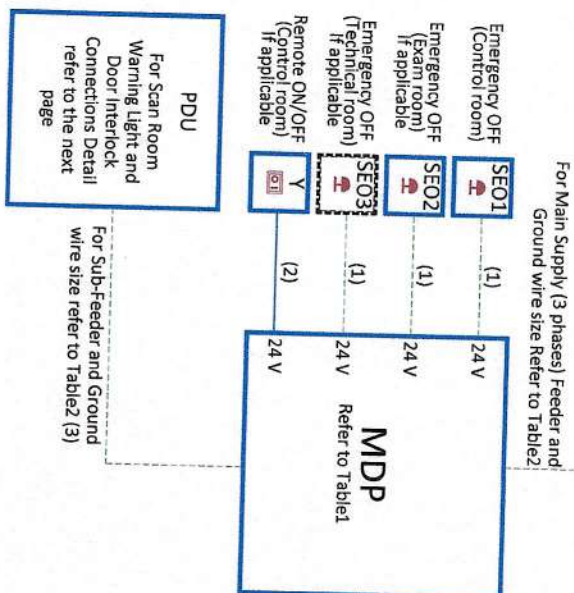


Table2:

Feeder Table

The information below assumes the use of copper wire, rated 75°C and run in steel conduit. All ampacity is determined in accordance with the National Electrical Code (NECA 70), Table 310-16 (2002). The ampacity of the circuit protection device listed above determines the minimum feeder size, except where total source regulation limits require a larger size. If the wire size does not match the above lists, please select the nearest wire size as per to local standards

Feeder length from Power Substation to MDP - ft (m)	Minimum Wire Size, AWG or MCM (mm²)/VAC							
	200 VAC	220 VAC	240 VAC	380 VAC	400 VAC	420 VAC	440 VAC	480 VAC
50 (15)	1/0 (55)	1/0 (55)	1/0 (55)	2 (35)	2 (35)	3 (30)	3 (30)	3 (30)
100 (30)	2/0 (70)	1/0 (55)	1/0 (55)	2 (35)	2 (35)	3 (30)	3 (30)	3 (30)
150 (46)	4/0 (100)	3/0 (85)	3/0 (85)	2 (35)	2 (35)	3 (30)	3 (30)	3 (30)
200 (61)	5/0 (125)	4/0 (100)	4/0 (100)	2 (35)	2 (35)	3 (30)	3 (30)	3 (30)
250 (76)	6/0 (170)	5/0 (125)	5/0 (125)	1 (45)	1 (45)	2 (35)	2 (35)	2 (35)
300 (91)	7/0 (215)	6/0 (170)	6/0 (170)	1 (45)	1 (45)	2 (35)	2 (35)	2 (35)
350 (107)	8/0 (275)	7/0 (215)	7/0 (215)	2/0 (70)	2/0 (70)	1/0 (55)	1/0 (55)	1/0 (55)
400 (122)	8/0 (275)	7/0 (215)	7/0 (215)	2/0 (70)	2/0 (70)	1/0 (55)	1/0 (55)	1/0 (55)
Sub-Feeder length from MDP to PDU - ft (m)								
32 (9.7536)	1 (45)	1/0 (55)	1/0 (55)	2 (35)	2 (35)	3 (30)	3 (30)	3 (30)

Run a dedicated 1/0 [50 mm²] or larger insulated copper ground wire from the power source to the MDP and from MDP to the PDU. Run the ground wire in the same raceway with the three-phase wires.

Notes:


- Wire size: 2x2mm² [14AWG] and 1x2mm² [14AWG] GND
- Power cable: 3 Meter/10', multi-conductor, 24V DC
- GE supplied MDP option E45021BB includes a 10 meter long power cable (H07RN-F) with wire size 4x50mm² and a 50 meter long control cable with wire size 2x1.5mm²

Table1:

GE Supplied Main Disconnect Panel (MDP)			
Region	CAT number	Amps	
Global except EMEA(440~480 V)	E4502BB	90	
Global except EMEA(380~420 V)	E4502BC	110	
EMEA(380~420 V)	E45021BB (3)	125	

--- Cable SUPPLIED BY CUSTOMER	
--- Cable SUPPLIED BY GE	
--- Equipment SUPPLIED BY CUSTOMER	
--- Equipment SUPPLIED BY GE	
MDP Main Disconnect Panel	
PDU Power Distribution Unit	
SEO Emergency OFF Button latching twist-to release style with two normally closed (NC) contacts	
Y System Remote Control with "ON" and "OFF" buttons	

SCAN ROOM WARNING LIGHT AND DOOR INTERLOCK

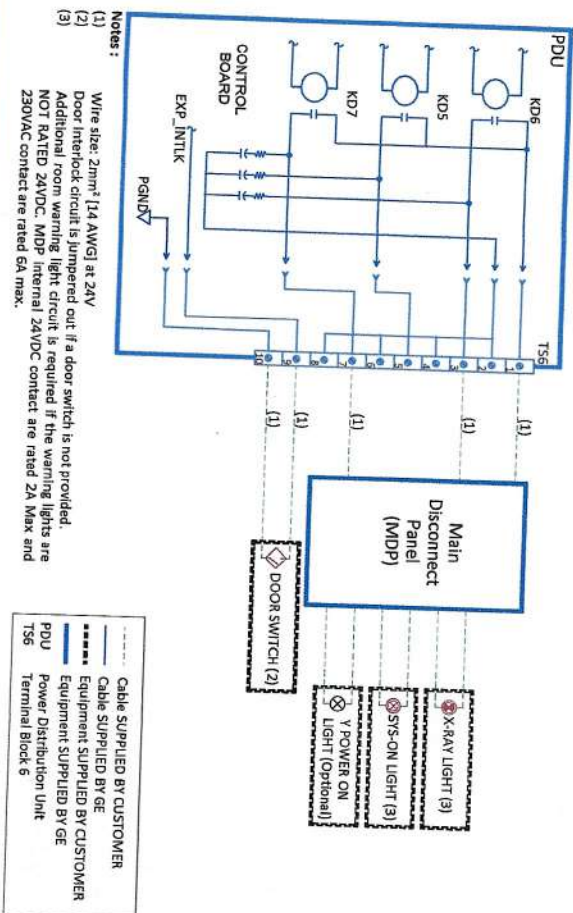

ENVIRONMENT

ALTITUDE

- The system shall meet all functional and performance specifications when placed in a room that is at an elevation of -150 m to 2,400 m [-492 ft to 7,875 ft] above sea level. For different altitudes refer to the PIM.

MAGNETIC FIELD SPECIFICATIONS

- Limit the magnetic interference to guarantee specified imaging performance.

GANTRY:

- Ambient static magnetic fields less than 1 Gauss.
- Ambient AC magnetic fields less than 0.01 Gauss peak.


OPERATOR CONSOLE:

- Ambient static magnetic fields less than 10 Gauss.

SYSTEM COMPONENT NOISE LEVEL

- Maximum Gantry Audible Noise Level:** The maximum ambient noise level is produced by the gantry during a CT scan acquisition. It is less than 70 dBA when measured at a distance of 1 m [3.3 ft] from the nearest gantry surface, in any direction.
- Maximum Console Audible Noise Level:** The maximum ambient noise levels is less than or equal to 54 dBA when measured 1 m [3.3 ft] up and 1 m [3.3 ft] away from the console at an ambient temperature of 26°C [79°F].

INTERCONNECTIONS

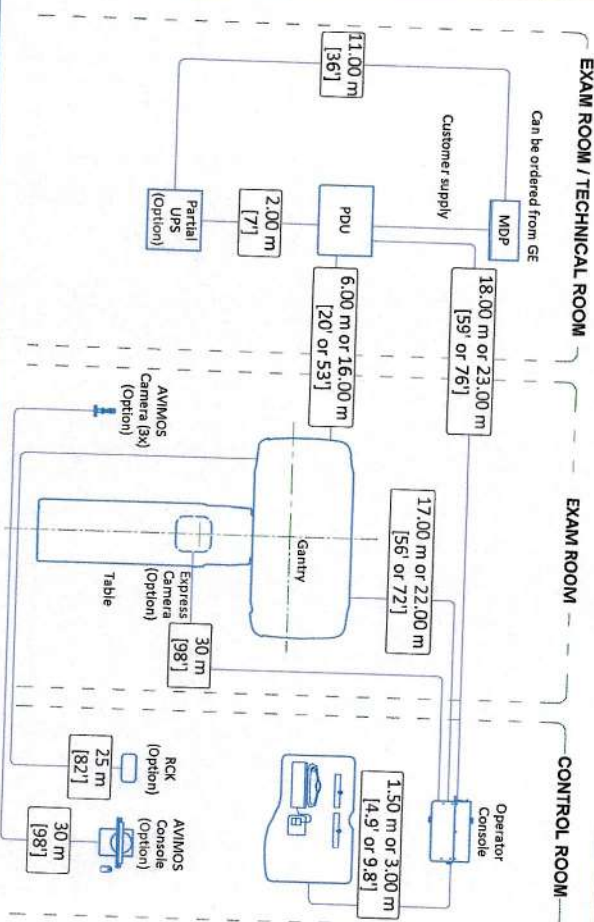

CONNECTIVITY REQUIREMENTS

Your new GE Healthcare imaging modality will require local and remote connectivity to enable our full range of digital support:


- Local connectivity** - This allows your system to connect to local devices such as PACS and modality worklist. We will require network information to configure the system(s), and a live ethernet port(s) prior to the delivery of the system(s).
- Remote connectivity** - Your GE Healthcare service warranty includes InSite™ (applicable to InSite capable products), a powerful broadband-based service which enables digital tools that can help guard your hospital against equipment downtime and revenue loss by quickly connecting you to a GE Healthcare expert.

Depending on product family and software version, imaging systems can be connected in one of the following methods:


1. TIS over TCP Port 443 (Preferred method for new products) via:
 - a. DNS resolution
 - b. Customer-provided Proxy or
 - c. GE Proxy (Available in some regions)
2. Site-to-Site IPsec VPN tunnel

Please provide the GE project manager with the contact information for the resource that can provide information required to set up these connections. GEHC will send out communication to these contacts, which will include the project's Connectivity requirements, and a Connectivity form. This form will need to be completed and returned to GEHC prior to delivery of the system to ensure the system is tested and connectivity is enabled prior to the completion of the installation.


GANTRY WITH GT1700V TABLE

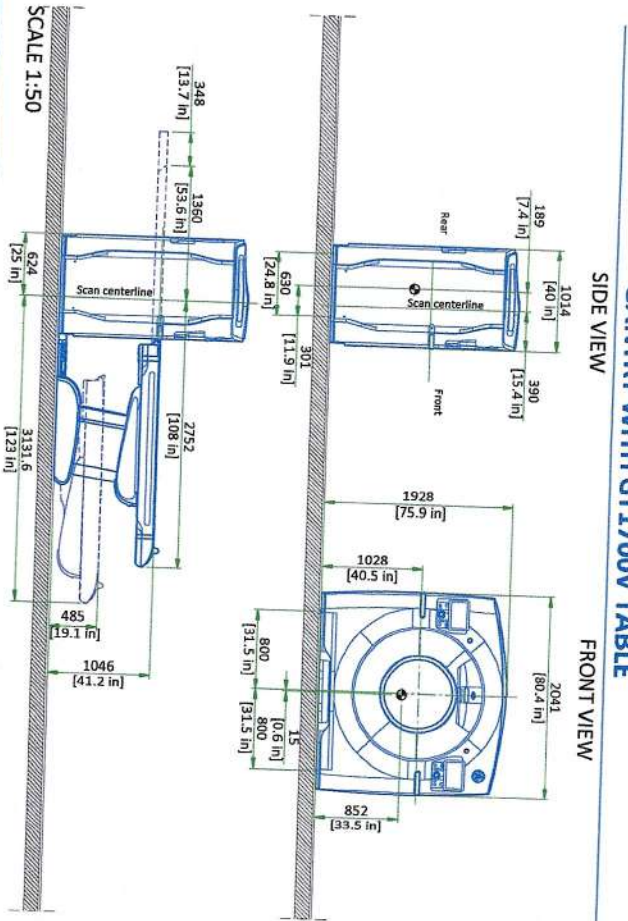


AURORA SWS TABLE

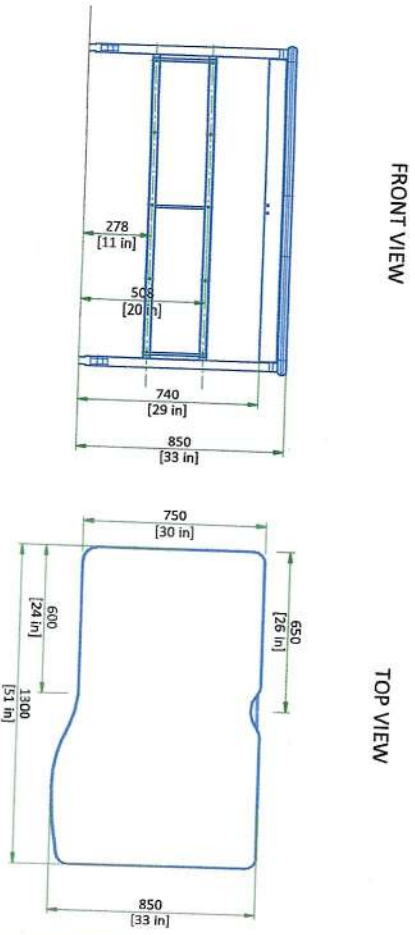

POWER DISTRIBUTION UNIT (PDU)

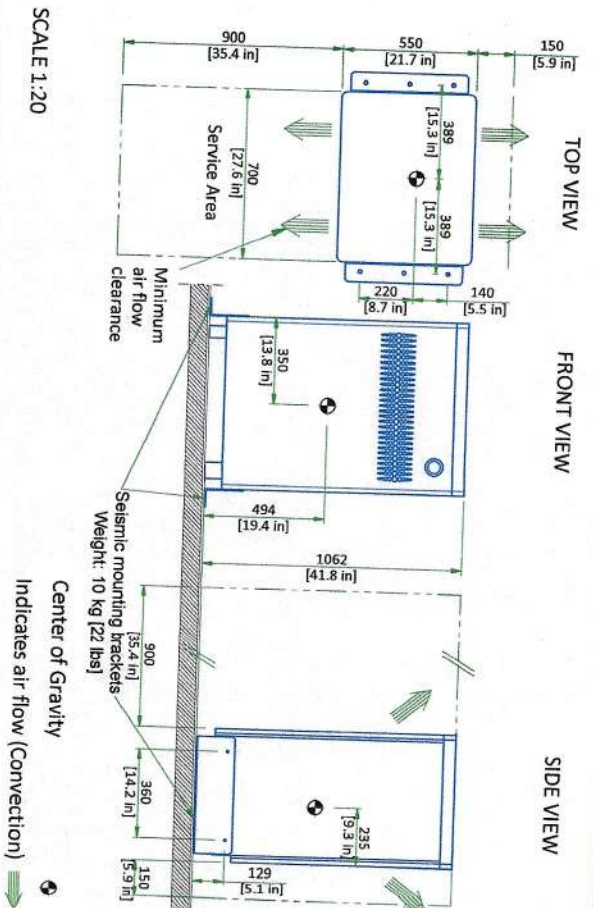

OPEN CONSOLE WITH Z8G4 HOST PC

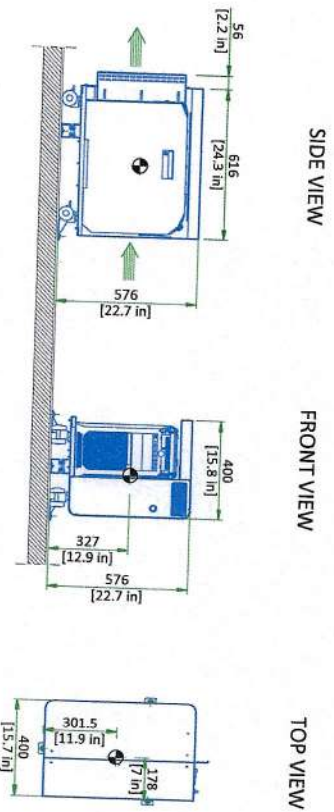

Weight: 40 kg [88 lbs]
SCALE 1:20

SUDABELT MEDICAL CO LTD

REVOLUTION MAXIMA

CT-8329269-FIN-00-A-DWG

Rev A/Date 24/FEB/2023

Equipment Dimensions (1)

09/12


Weight: 64.2 kg [142 lbs]
SCALE 1:20

Center of Gravity
Indicates air flow (Convection)


CAUTION
Mechanical Engineering Ltd.

MAIN DISCONNECT PANEL


TOP VIEW

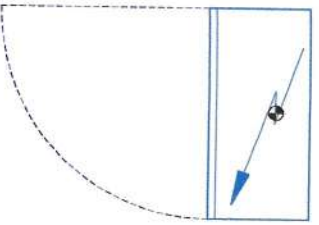

FRONT VIEW

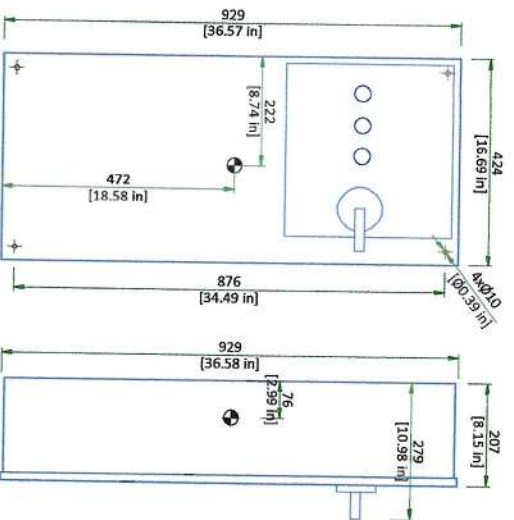

SIDE VIEW

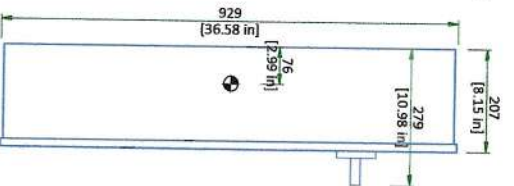

Center of gravity ⬤
SCALE 1:10

Handwritten signature and date: 24/02/2023

TEMPERATURE AND HUMIDITY SPECIFICATIONS

IN-USE CONDITIONS

	EXAM ROOM			CONTROL ROOM			TECHNICAL ROOM		
	Min	Recommended	Max	Min	Recommended	Max	Min	Recommended	Max
Temperature	18°C	22°C	26°C	18°C	22°C	26°C	18°C	22°C	26°C
	64°F	72°F	79°F	64°F	72°F	79°F	64°F	72°F	79°F
Relative humidity (1)	30% to 60%			30% to 60%			30% to 60%		

STORAGE CONDITIONS

Temperature	0°C to +30°C 32°F to +86°F
Temperature gradient	≤ 3°C/h
Relative humidity (1)	≤ 51°F/h up to 70%
Humidity gradient	≤ 5%/h

Storage longer than 6 months is not recommended.
(1) Non-condensing

AIR RENEWAL

According to local standards.

NOTE
In case of using air conditioning systems that have a risk of water leakage it is recommended not to install it above electric equipment or to take measures to protect the equipment from dropping water.

HEAT DISSIPATION DETAILS

ROOM	DESCRIPTION	Max (kW)	Max (BTU)
Exam Room	Gantry	5.48	18700
	GT1700V/GT2000 Table	0.3	1030
	TOTAL	5.78	19730
Exam Room or Technical Room*	Power Distribution Unit	1.0	3400
	TOTAL	1.00	3400
Control Room	Operator console	0.84	2860
	LCD monitor (Total amount of 2 monitors)	0.1	340
	TOTAL	0.94	3200

*Technical Room is not mandatory, the placements of those elements are recommended in the Exam Room.

DELIVERY


THE CUSTOMER/CONTRACTOR SHOULD:

- Provide an area adjacent to the installation site for delivery and unloading of the GE equipment.
- Ensure that the dimensions of all doors, corridors, ceiling heights are sufficient to accommodate the movement of GE equipment from the delivery area into the definitive installation room.
- Ensure that access routes for equipment will accommodate the weights of the equipment and any transportation, lifting and rigging equipment.
- Ensure that all necessary arrangements for stopping and unloading on public or private property not belonging to the customer have been made.

DIMENSIONS OF DELIVERY WITH DOLLY TRANSPORT EQUIPMENT

EQUIPMENT	DIMENSIONS			WEIGHT	
	LENGTH	WIDTH	HEIGHT		
GANTRY	2775 mm	109 in	1291 mm	1835 kg	4049 lbs
			2000 mm		
			2489 mm		
GT1700V TABLE			762 mm		
			1143 mm	576 kg	1270 lbs
			45 in		

GANTRY DELIVERY

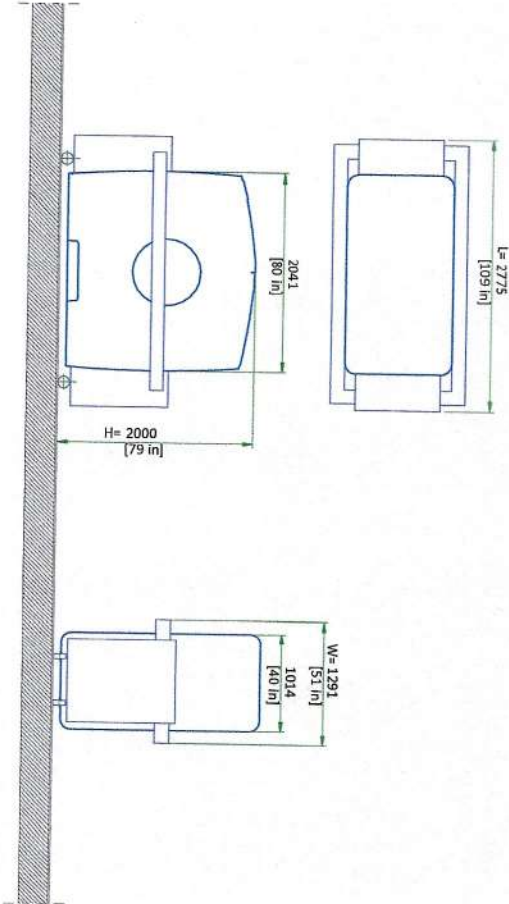

- The gantry is shipped on a dolly equipped with elevating casters (normal shipping configuration).
NOT TO SCALE

DISCLAIMER

GENERAL SPECIFICATIONS


- GE is not responsible for the installation of developers and associated equipment, lighting, cassette trays and protective screens or derivatives not mentioned in the order.
- The final study contains recommendations for the location of GE equipment and associated devices, electrical wiring and room arrangements. When preparing the study, every effort has been made to consider every aspect of the actual equipment expected to be installed.
- The layout of the equipment offered by GE, the dimensions given for the premises, the details provided for the pre-installation work and electrical power supply are given according to the information noted during on-site study and the wishes expressed by the customer.
- The room dimensions used to create the equipment layout may originate from a previous layout and may not be accurate as they may not have been verified on site. GE cannot take any responsibility for errors due to lack of information.
- Dimensions apply to finished surfaces of the room.
- Actual configuration may differ from options presented in some typical views or tables.
- If this set of final drawings has been approved by the customer, any subsequent modification of the site must be subject to further investigation by GE about the feasibility of installing the equipment. Any reservations must be noted.
- The equipment layout indicates the placement and interconnection of the indicated equipment components. There may be local requirements that could impact the placement of these components. It remains the customer's responsibility to ensure that the site and final equipment placement complies with all applicable local requirements.
- All work required to install GE equipment must be carried out in compliance with the building regulations and the safety standards of legal force in the country concerned.
- These drawings are not to be used for actual construction purposes. The company cannot take responsibility for any damage resulting therefrom.

CUSTOMER RESPONSIBILITIES

- It is the responsibility of the customer to prepare the site in accordance with the specifications stated in the final study. A detailed site readiness checklist is provided by GE. It is the responsibility of the customer to ensure all requirements are fulfilled and that the site conforms to all specifications defined in the checklist and final study. The GE Project Manager of Installation (PMI) will work in cooperation with the customer to follow up and ensure that actions in the checklist are complete, and if necessary, will aid in the rescheduling of the delivery and installation date.
- Prior to installation, a structural engineer of record must ensure that the floor and ceiling is designed in such a way that the loads of the installed system can be securely borne and transferred. The layout of additional structural elements, dimensioning and the selection of appropriate installation methods are the sole responsibility of the structural engineer. Execution of load bearing structures supporting equipment on the ceiling, floor or walls are the customer's responsibility.

RADIO-PROTECTION

- Suitable radiological protection must be determined by a qualified radiological physicist in conformation with local regulations. GE does not take responsibility for the specification or provision of radio-protection.


THE UNDERSIGNED, HEREBY CERTIFIES THAT I HAVE READ AND APPROVED THE PLANS IN THIS DOCUMENT.		
DATE	NAME	SIGNATURE
18/10/2023	Ogo Kayode	

CUSTOMER SITE READINESS REQUIREMENTS

REQUIRED MANUALS FOR SYSTEM PRE-INSTALLATION

Description	Document Number*
Product specific Pre-Installation Manual	Refer to cover page
*documents can be accessed in multiple languages at https://customer-docs.cloud.gelifehealthcare.com/#/cdp/dashboard	

- A mandatory component of this drawing set is the GE Healthcare Pre-Installation manual. Failure to reference the Pre-Installation manual will result in incomplete documentation required for site design and preparation.
- The items on the GE Healthcare Site Readiness Checklist **DOC1809666** are REQUIRED to facilitate equipment delivery to the site. Equipment will not be delivered if these requirements are not satisfied.
- Any deviation from these drawings must be communicated in writing to and reviewed by your local GE Healthcare installation project manager prior to making changes.
- Make arrangements for any rigging, special handling, or facility modifications that must be made to deliver the equipment to the installation site. If desired, your local GE Healthcare installation project manager can supply a reference list of rigging contractors.
- New construction requires the following:
 1. Secure area for equipment,
 2. Power for drills and other test equipment,
 3. Restrooms.
- Provide for refuse removal and disposal (e.g. crates, cartons, packing)
- For CT systems it is required to minimize vibrations within the scan room. It is the customer's responsibility to contract a vibration consultant/engineer to implement site design modifications to meet the GE vibration specification. Refer to the system Pre-Installation manual for vibration specifications.

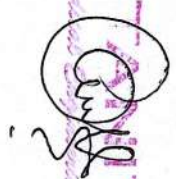
Ogo Kayode