


T-2100 treadmill

Advanced performance in exercise testing

Advances incorporated into the T-2100 treadmill place it at the forefront of patient safety, dependability and overall performance – the culmination of continual technological development and attention to detail specially designed to meet the high quality standards you expect from GE.

The T-2100 treadmill's extra-long, 152 cm (60-inch) treadmill surface provides an extra margin of safety for patients. Its zero start feature slowly builds up speed for safe and smooth operation. And it can be slowed to a gradual stop or stopped immediately with the standard emergency stop button. Full front and side handrails are also standard on the T-2100 treadmill.

With fewer moving parts than other treadmills, the T-2100 treadmill delivers extremely smooth, quiet operation that facilitates accurate blood pressure measurements even at high workloads. Thanks to digital microprocessor control, its settings never need to be calibrated. And a 204 kg (450-pound) weight capacity and speed range of 0-21,7 km/h (0-13.5 miles per hour) set the Series 2000 apart from other treadmills.

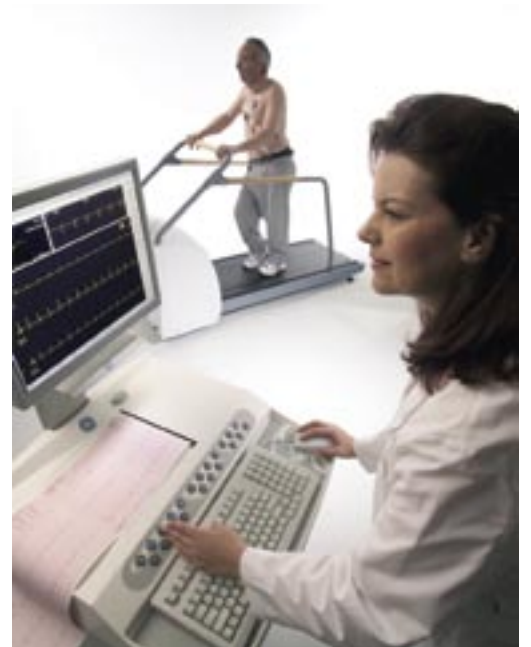

Feature highlights

Patient safety

- 152 cm (60-inch) walking surface
- Two stopping modes
- Standard full handrail set
- Standard emergency stop switch
- Starts from 0 km/h (mph)
- 0 to 21,7 km/h (0 to 13.5 mph) speed range
- Weight capacity of 204 kg (450 lbs.)
- 0 to 25 percent grade range

Reliability

- Enhanced tolerance to electrical changes for greater patient safety
- 3.0 horsepower motor
- Self-calibrating for speed and grade
- Digital control insures precise speed and grade control

Interface

- CASE® Cardiac Assessment System for Exercise
- MAC 1200 ST, MAC 5000 and MAC 5500 ECG acquisition system
- CardioSoft multifunction PC-based diagnostic software
- Also interfaces with many legacy products

T-2100 treadmill	
Physical specifications	
Weight	181.4 kg (400 lbs.)
Walking area	45.7 x 152.4 cm (18 x 60 in.)
Floor space required	73.7 x 195.6 cm (29 x 77 in.)
Walking surface height	14.0 cm (5.5 in.)
Handrail height above walking surface	Front – 101.6 cm (40 in.) Side – Maximum 88.9 cm (35 in.) descending at a 5° angle to 81.3 cm (32 in.)
Environmental specifications	
Power requirements	200 to 240 VAC, 50/60 Hz, single-phase, 20 amperes, NEMA 6-20 R wall socket (dedicated circuit recommended)
Operating conditions	
Ambient temperature	10° C to 40° C (50° F to 95° F)
Relative humidity	15% to 90%
Storage conditions	
Temperature range	-18° C to 49° C (0° F to 120° F)
Relative humidity	15% to 90%
Performance specifications	
Maximum rated load	204 kg (450 lbs.)
Belt speed range	0 to 21.7 km/h (0 to 13.5 mph) at 60 Hz continuously variable; 0 mph start-up
Driver motor	3.0 hp
Elevation motor	1/10th hp
Elevation range	0 to 25% grade, continuous variable
Maximum leakage current	100µA, chassis to ground

GE Medical Systems Information Technologies GmbH
Munzinger Strasse 3
79111 Freiburg, Germany
Tel. +49 761 4543 0 • Fax +49 761 4543 233

www.gehealthcare.com

©2005 General Electric Company – All rights reserved.

General Electric Company reserves the right to make changes in specifications and features shown herein, or discontinue the product described at any time without notice or obligation.

GE, GE Monogram and CASE® are trademarks of General Electric Company.

GE Medical Systems Information Technologies GmbH, a General Electric company, going to market as GE Healthcare.



GE imagination at work